

Fairworth Chiropractic Tables



Fairworth Electric Chiropractic Table Series

Main features

The Fairworth chiropractic table is a very practical manual chiropractic table with flexion, traction and drop functions. In particular, its vertical and lateral lumbar flexion function is very practical for chiropractors. It is equipped with basic head, thoracic, lumbar and pelvic drop mechanisms, and the cervical and lumbar traction guide rail. Its heavy-duty frame, Stand-by Castor System, Walk-Free Control System, Self-Skin Injection Foam and Rotating Foot Support Extender all mark Fairworth's ability to become the best chiropractic table in the world.

The great Fairworth chiropractic table---the most versatile chiropractic and treatment table in the world;



Video




Video



Self-Skin Injection Foam

Our self-skinning sponge pad is made of 200kg high-density polyether through high-pressure sealed injection molding, and the skin is integrally formed. It is featuring durable, water and oil resistance.

Master of Chiropractic Setting

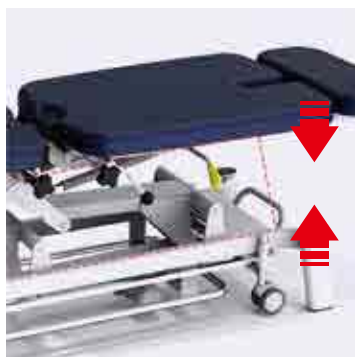

Powerlift Height Adjustment

Adjust your table to the ideal height with the tap of a footswitch. The table has an elevation range of 55-97cm/21.6"-38", is rated to 200 kg. (440 lbs.)



Flippable Headpiece

The headpiece upper layer can flip up 0°-30° in reverse direction, supported by pneumatic pump.



Tilting Adjustable Headpiece

The Fairworth features a -40°~30° comfortable and easily pneumatic adjustable tilting headpiece.



Straight Cervical Drop with Adjustable Tension

Our headpiece features a crisp and clean drop mechanism and an adjustable tension knob.



Manual Headpiece Traction

The Headpiece traction perfectly mobilize and tract the neck, cervical and spine to reach the goal of full relaxation.



Height Adjustable Front Arms

The height adjustable front armrests can adjust the height in the range of 0-20 cm and can be locked at any position.



Straight Thoracic Drop with Adjustable Tension

The Fairworth features a thoracic drop mechanism and an adjustable tension knob.




Straight Lumbar Drop with Adjustable Tension

The Fairworth features a lumbar drop mechanism and an adjustable tension knob.




Tension Adjustable Pelvic Drop

The Fairworth features a pelvic drop mechanism and an adjustable tension knob.

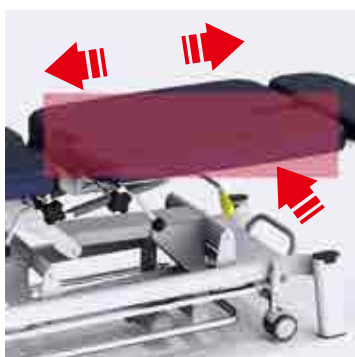

Manual Lumbar Traction

Manual lumbar traction mechanism with built-in guide rail provides smooth traction feeling.

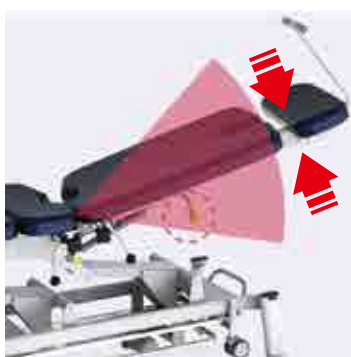

Manual Lumbar Vertical Flexion

The smooth vertical flexion range is -40°~20°. The pneumatic supporting pump can be locked at any position while flexed.




Manual Lumbar Lateral Flexion

The lateral flexion range is 18° to either side of midline. It can be locked when not needed.

Master of Chiropractic Setting



Flippable Foot Plate

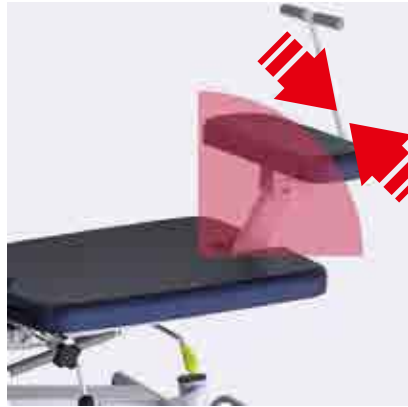

Rotating Foot Support Extender

Foot Support Extender adds an additional 8 inches. It can be used with Flippable Foot Plate to fix the ankle of patients.



Removable Flexion Control Handle

The removable flexion control handle will be equipped on Fairworth-360 and Fairworth-380 which have the flexion mechanism.

Meaning of Product Code

Fairworth product code consists of three model codes, four option codes. The four option codes for Fairworth table depend on the configuration of various options available. Each position in the code has a specific meaning, and corresponds to one or more characteristics of the table. Our standard model of Fairworth itself does not contain Walk-free Control System, Stand-by Castor System, Electrical Safety Lock and body Straps, all optional parts need to be purchased separately. All standard models adopt front armrests, powerlift height adjustment, self-skin sponge system.

POSITION CODE 1~3

Product Standard Model

POSITION CODE 4

Control System

- A - Single Pedal
- B - Hand Control
- C - Walk-free Controller
- D - Walk-free Controller + Hand Control

POSITION CODE 5

Stand-by Castor System

- 1 - With Stand-by Castor System
- 0 - Without Stand-by Castor System

POSITION CODE 6

Safety Lock System

- 1 - With Electrical Safety Lock
- 0 - Without Electrical Lock System

POSITION CODE 7

Body Straps

- 1 - With 2 Body Straps
- 0 - Without Body Strap

Model: 3 8 0 - C 1 1 1

Code Position: 1 2 3 - 4 5 6 7

Model Codes

Option Codes



Stand-by
Castor System

4 Control System

6 Safety Lock System

7 Body Straps



POSITION 4

			-				
1	2	3		4	5	6	7

Control System

We provide a variety of motor control options, all controllers must be installed in advance by original factory. Please indicate when ordering.



A Single Pedal

Standard Controller of Fairworth table. It is used to raise or lower the height of the table.



B Hand Control

2 Button Hand Control is used for adjusting the height of the table. It is applicable for all Fairworth models.



C Walk-Free Control System

Walk-Free Control makes height adjustment easily from all the way around the bed. It is applicable to all Fairworth Models.



D Walk-Free Control + Hand Control

Walk-Free Control is for adjusting the table's height. 2 Button Hand Control is for raising/lowering the height of the table. It is applicable for all models.

POSITION 5

			-				
1	2	3		4	5	6	7

Stand-by Castor System



The lift able Stand-by Castor System adopts the medical castors and the connecting rod operating system. The four castors can be released and retracted synchronously by foot pedal control. The bed can be moved easily when wheels are released, and the stability of the bed will be increased when the castors are folded up. Castor system must be installed by the original factory. Users are not allowed to install the castor system without permission to avoid accidental injury.

1 With Stand-by Castor System

0 Without Stand-by Castor System

POSITION 6

			-				
1	2	3		4	5	6	7

Safety Lock System

The safety lock adopts the traditional mechanical lock, which is safe and reliable. It is used to completely cut off the power supply of the motor and effectively avoid the operation and accidents of irrelevant personnel.




1 With Electrical Safety Lock

According to the European law, the electric lifting bed must be equipped with safety devices to avoid accidental injury caused by the operation of children or non-designated personnel. Each lockset is equipped with two keys.

0 Without Electrical Lock System

POSITION 7

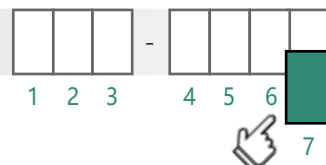

Body Straps




1 With 2 x Body Straps

To provide better restraint for patients, fixed by velcro, there are guide rails on both sides of the table frame, which can fix the straps at the required position.

0 Without Electrical Lock System

Standard Model

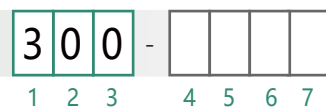

Fairworth-300

Fairworth-300 is a powerlift height adjustable bed with straight lumbar drop, pelvic drop, headpiece drop, thoracic drop mechanisms and all drop mechanisms are tension adjustable by adjusting knobs. The standard model includes the front arms, foot support extender and tilting headpiece. The manual traction, the vertical and lateral flexion mechanisms and flexion control handle are not included.



Video

Specification:

- Working Weight: 440lbs/200kg
- Table Length: 69"/175cm
- Width: 22"/56cm
- Height range: 21.6"/55cm~38"/97cm
- Weight: 252lbs/115kg

Standard Features:

- 5 Year Limited Warranty
- Heavy duty metal frame
- Self-skining injection foam table top
- Powerlift table height adjustable
- Tilting adjustable and flappable Headpiece.
- Headpiece straight drop with adjustable tension
- Thoracic straight drop with Adjustable Tension
- Pelvic straightdrop with Adjustable Tension
- Lumbarstraight drop with Adjustable Tension
- Rotating Foot Support Extender

Standard Model

3 6 0 -

1 2 3 4 5 6 7

Fairworth-360

Fairworth-360 is a powerlift height adjustable flexion and traction table, equipped with manual lumbar vertical and lateral flexion mechanism, manual lumbar traction mechanism and manual head traction mechanism. The standard model includes the removable flexion control handle, front arms, foot support extender and tilting headpiece. The lumbar drop, pelvic drop, headpiece drop, thoracic drop are not included.



Standard Features:

- Manual lumbar vertical flexion
- Manual lumbar lateral flexion
- Manual lumbar traction
- Manual headpiece traction
- Rotating Foot Support Extender
- Heavy duty metal frame
- Self-skining injection foam table top
- Powerlift table height adjustable
- Tilting adjustable and flappable Headpiece

Specification:

- Working Weight: 440lbs/200kg
- Table Length: 69"/175cm
- Width: 22"/56cm
- Height range: 21.6"/55cm~38"/97cm
- Weight: 252lbs/115kg

Standard Model

3 8 0 -

1 2 3 4 5 6 7

Fairworth-380

Fairworth-380 is a full-function powerlift height adjustable flexion, traction and drop table, equipped with manual lumbar vertical and lateral flexion mechanism, manual lumbar traction mechanism, manual head traction mechanism and straight lumbar, pelvic, headpiece, thoracic drop mechanisms. The standard model includes the removable flexion control handle, front arms, foot support extender and tilting headpiece.



Specification:

- Working Weight: 440lbs/200kg
- Table Length: 69"/175cm
- Width: 22"/56cm
- Height range: 21.6"/55cm~38"/97cm
- Weight: 285lbs/130kg

Standard Features:




- 5 Year Limited Warranty
- Heavy duty metal frame
- Self-skining injection foam table top
- Powerlift table height adjustable
- Tilting adjustable and flappable Headpiece.
- Manual lumbar vertical flexion
- Manual lumbar lateral flexion
- Manual lumbar traction
- Manual headpiece traction
- Headpiece straight drop with adjustable tension
- Thoracic straight drop with Adjustable Tension
- Pelvic straight drop with Adjustable Tension
- Lumbar straight drop with Adjustable Tension
- Rotating Foot Support Extender

Customize Your Own Table

			-				
1	2	3		4	5	6	7

Determine Product Code

☒ Basic Configuration
 ☐ Optional Configuration
 ☒ Unavailable

Product code	300	360	380
Image			
- Length	68.89"/175cm		
- Width	22.04"/56cm		
- Weight	286.6lbs/130kg		
- Height	21.65"/55cm~38.18"/97cm		
- Working Load	440.92lbs/200kg		
Powerlift Table Height Adjustment	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Flippable Pneumatic Headpiece	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Pneumatic Tilting Headrest	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Manual Lumbar Vertical Flexion	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Manual Lumbar Lateral Flexion	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Lumbar Traction	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Headrest Traction	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Removable Flexion Control Handle	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Straight Cervical Drop	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Straight Thoracic Drop	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Straight Lumbar Drop	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Straight Pelvic Drop	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Position Code 4 Control System			
A - Single Pedal	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
B - Hand Control	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
C - Walk free control	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
Position Code 5 Standby Castor System			
1 - With Standby Castor System	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Position Code 6 Safety Lock System			
1 - With Safety Lock System	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Position Code 7 Body Straps			
1 - With Body Straps	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

Option Varieties

Fairworth-300



300-A010



300-A110



300-C110

Fairworth-360



360-A011



360-A111



360-C111

Fairworth-380



380-A011



380-A111



380-C111



Biocure Single Member P.C.
Epidavrou 2, Evosmos, 56224 Thessaloniki, Greece